







Distillation Unit K-350 light

Powerful steam distillation for versatile applications




- Ease of use
Straightforward operation due to intuitive programming
- Best accuracy results
The slim glass splash protector reduces dead volume and facilitates vapour transfer
- Efficient and proven
The powerful steam generator allows fast analyte separation from matrix
- Versatile applications
Easy implementation of steam distillation into your workflow

	Matrix	Analyte	Norm / Regulation
	Alcoholic beverages (e.g. beer, wine, spirit, liqueur)	Alcohol	EEC 2870/2000 (spirits) EEC 2676/1990 (wine) OIV-MA-AS312-01A (wine)
	Alcoholic and fruit-based beverages (e.g. cider, juice)	Volatile acids	EEC 2676/1990 (wine)
	Leather	Chlorophenol	DIN EN ISO 17070
	Plants	Essential oils	–

Distillation Unit K-350 light Features



Powerful distillation
Proven steam generator for
numerous applications

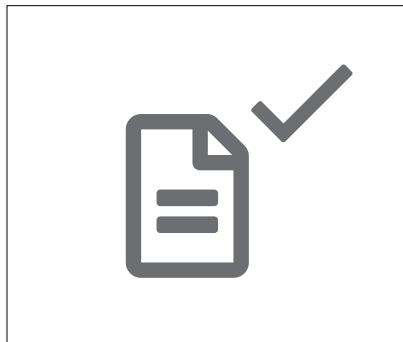

Versatility
Multiple combinations with
analysers such as titrators and
densimeters




Save tap water
Start and stop cooling water
during the distillation process or
use recirculating chiller



Safety features
Protective sensors for sample
tube, front and service door
ensure safe operation



Accessories
Choose the appropriate accessory
for your application

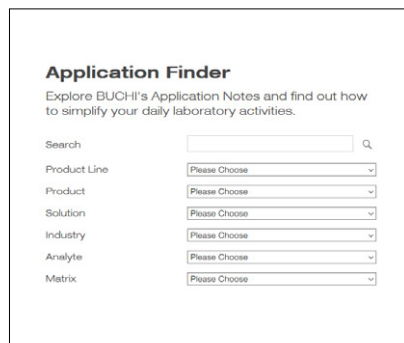

Application finder
Find your application for a quick
start

“The introduction of the BUCHI Distillation Unit K-350 in our company made a fast and easy alcohol determination possible. It makes the alcohol determination much more reliable and not dependent from the executor anymore.”

Dr. Patrick van Zuidam, Managing Director Zuidam Distillers BV, Netherlands

Quality in your hands

BÜCHI Labortechnik AG
CH - 9230 Flawil
T +41 71 394 63 63
F +41 71 394 64 64
info@buchi.com

www.buchi.com

